



MADE IN GERMANY

Wassermann is a medium-sized, owner-operated company with head office in Hamburg. We design and produce high-quality technical devices in Germany. Our main users are dental laboratories, dental clinics and various industrial areas worldwide.

Why a Wassermann product?


Wassermann products are manufactured from high-quality materials and the proportion of plastics used is intentionally kept as low as possible. Our devices are known for their suitability for daily use, functionality and longevity. As a result, the users place high levels of trust in the Wassermann brand. We feel it is our duty to support our customers on a daily basis when it comes to safety, service and productivity. Technological perfection, optimum quality and especially sustainability are the benchmarks we set. Thanks to the combination of these values, the family-owned company, since its establishment in 1927, has gained worldwide recognition as a manufacturer of high-quality products.



Wapo-Ex 12 II



Boil Out Unit

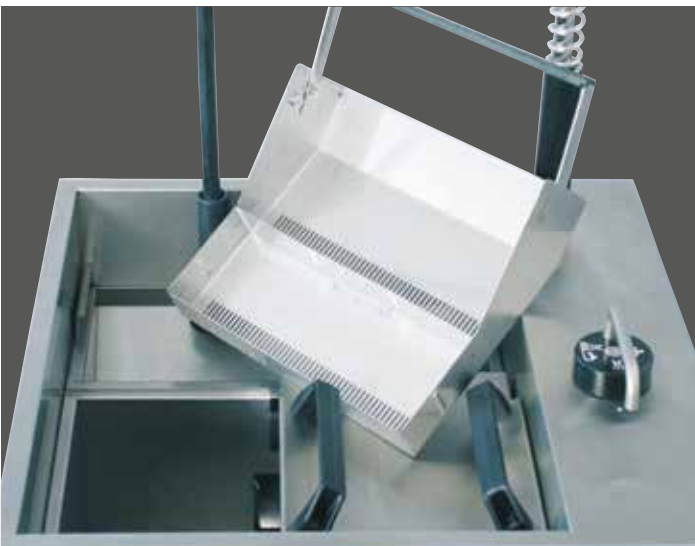

Further information are available on our website at www.wassermann.hamburg

Your Dental Dealer



805 Marathon Parkway Suite 130
Lawrenceville, GA 30046
www.zublerusa.com

GA Office (770) 921-2131
TX Office (972) 600-9187
sales@zublerusa.com



Boil Out Unit **Wapo-Ex 12 II**

The boil out unit **Wapo-Ex 12 II** stands out thanks to its efficient applications planning, programmable timer and automatic wax removal. This free standing unit with integrated spray gallows scalds up to 12 flask halves simultaneously and fully automatically. A touch screen offers diverse functions and reproducible procedures by simple programming. An elaborate boil out unit manufactured using high quality materials and components and designed for durability.

Integrated spray gallows
for manual scalding of small parts

Selector lever for different functions as
automatic scalding, manual scalding
(hand spray) and pump down

Touch screen for
fast and easy operation

Energy saving
thermal insulation

High quality construction,
most parts made of stainless steel



Major Features

- Free standing unit for scalding up to 12 flask halves
- Efficient usage thanks to its programmable timer and automatic wax removal
- Safety features with visual and audible signals, e.g. electronic water level gauge, leakage signal, inlet limitation
- Energy saving because of the intelligent heating control
- **Included in delivery:**
 - Spray gallows
 - Service set (1 L lime remover and 2 L wax remover)
 - Supply hose with water inlet valve and gasket
 - Discharge hose with hose band clamp Ø 12-20 x 9
 - Spray arm tube with seal
 - 2 flask baskets
 - 2 covers

Accessories

Service set (1 L lime remover and 2 L wax remover)

Item no.

170330

Wapo-Ex 12 II

Voltage

220-240 V / 50/60 Hz

Power consumption

12 A

Output

max. 3300 W

W x H (without spray gallows) x D

650 x 900 x 620 mm

W x H (with spray gallows) x D

650 x 1600 x 620 mm

Weight

87 kg

Sound level

≤ 70 dB (A)

Water capacity

approx. 50 l

Water feed rate

4.5 l / min

Water drain height

470 mm

