





Hand spray for manual scalding

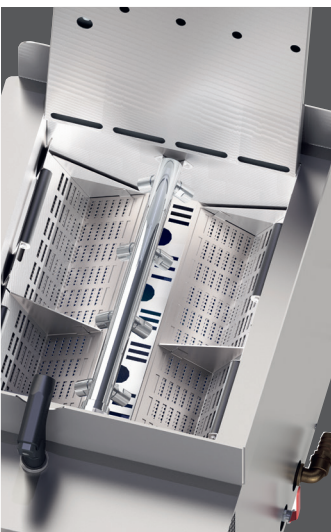

Convenient boil out process, easy access to flask baskets

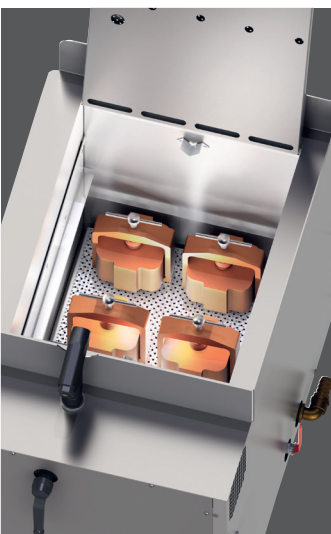

Reliable polymerisation results

### MADE IN GERMANY

Wassermann is a medium-sized, owner-operated company with head office in Hamburg. We design and produce high-quality technical devices in Germany. Our main users are dental laboratories, dental clinics and various industrial areas worldwide.

### Why a Wassermann product?

Wassermann products are manufactured from high-quality materials and the proportion of plastics used is intentionally kept as low as possible. Our devices are known for their suitability for daily use, functionality and longevity. As a result, the users place high levels of trust in the Wassermann brand. We feel it is our duty to support our customers on a daily basis when it comes to safety, service and productivity. Technological perfection, optimum quality and especially sustainability are the benchmarks we set. Thanks to the combination of these values, the family-owned company, since its establishment in 1927, has gained worldwide recognition as a manufacturer of high-quality products.



Further information are available on our website at [www.wassermann.hamburg](http://www.wassermann.hamburg)




/zublerusallc

Your Dental Dealer



[www.zublerusa.com](http://www.zublerusa.com)  
[sales@zublerusa.com](mailto:sales@zublerusa.com)

### Wapo-Ex 8 Plus

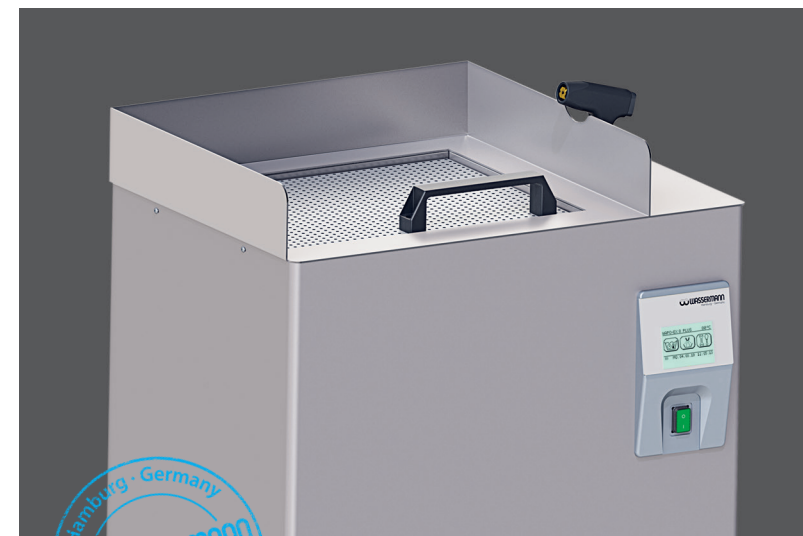

Boil Out and Polymerisation Unit



Wassermann Dental-Maschinen GmbH  
Rudorffweg 15-17 · 21031 Hamburg · Germany  
Phone: +49 (0)40 730 926-20/ -24 · Fax: +49 (0)40 730 37 24  
[export@wassermann-dental.com](mailto:export@wassermann-dental.com) · [www.wassermann.hamburg](http://www.wassermann.hamburg)  
[www.facebook.com/WassermannDentalmaschinen](https://www.facebook.com/WassermannDentalmaschinen)

# Combination Automat

## Wapo-Ex 8 Plus

### Boil Out and Polymerisation Unit

The **Wapo-Ex 8 Plus** enables the automatic boil out of up to 8 flask halves or the polymerisation of up to 8 flasks in a clamp. The relevant parameters are at once visible on the touch screen, the input can be made intuitively. The user can adjust the numerous parameters and options (such as week timer) to his individual requirements by simple programming.



### Major Features

- Compact tabletop automat for boil out up to 8 flask halves or polymerisation up to 8 flasks in a clamp
- Thermostat protected, gasketless submersible pump, maintenance-free
- Spray system in the middle, easy to mount/ dismount
- Quick exchange of both functions boil-out and polymerisation
- Optical and acoustic lack of water warning
- Easy to clean without tools
- **Included in delivery:**
  - Service set (4 filter mats, 1 L lime remover, 2 L wax remover)
  - Discharge hose, 1.5 m
  - Hose clip Ø 20-32 x 12.2

Accessories	Item no.
Mobile stand (W 620 x H 550 x D 430 mm)	170120
Service set (4 filter mats, 1 L lime remover, 2 L wax remover)	170021
Discharge hose per meter	320030

Wapo-Ex 8 Plus	Item no. 170977
Voltage	220–240 V / 50/60 Hz
Power consumption	9 A
Output	2000 W
W x H x D max. W x H x D (housing)	660 x 585 x 480 mm 615 x 550 x 425 mm
Weight	41 kg
Sound level	≤ 70 dB (A)
Water capacity	21 l