



MPF

BRUSH CO



Many people say we have a think outside the box aproach,
we like to think there is no box at all.

GEORGE PAPPS
CEO

At MPF Brush Co. our primary goal is to continue to improve our customer's performance by offering a unique range of laboratory hand tools, brushes, pal-ettes, accessories and rotary instruments at a reasona-ble price, without compromising on quality.

With more than 20 years in the market we are the leaders in premium dental instruments worldwide.

Collaborating with opinion leading dental technicians from around the world has enabled us to 'think outside the box' when designing unique dental ceramic hand tools.

Thanks to our loyal customers, MPF Brush Co. has been producing dental instruments for some time now and our product range is rapidly growing. It is you that drives us to create new and innovative tools that are not only beautiful and ergonomic in design, but also functional and affordable.

We'd like to take this opportunity to thank all of our customers worldwide for their ongoing support.

01 PG 6 - 11**MPF Hybrid Series**

MPF Hybrid Brush Tips	(8-9)
MPF Hybrid Plus Series	(10-11)

02 PG 12-35**Brushes & Brush Kits**

Optimum™ Mens Brushes & Brush Kits	(14)
Optimum™ Ladies Brushes & Brush Kits	(15)
Optimum Master Brushes & Brush Kits	(17)
MPF Revolution Brushes & Brush Kits	(18)
MPF Revolution Master Brushes & Brush Kits	(19)
MPF Revolution Handles	(21)
Jungo Endo Signature Brushes	(22)
MPF Signature Brush	(22)
Stipple by Stefan Picha	(22)
MPF Stadar Brushes & Brush Kits	(23)
MPF Zirconia Stain Brushes & Brush Kit	(24)
MPF Synthesis Brushes & Brush Kit	(25)
MPF Classic Kolinsky Brushes & Brush Kit	(27)
Replacement Tips	(32-33)

Dental Surgical

MPF Composite Stain Brushes & Brush Kit (Straight)	(28)
MPF Composite Stain Brushes & Brush Kit (Angled)	(28)
MPF Composite Flat Brushes & Brush Kit (Straight)	(29)
MPF Composite Flat Brushes & Brush Kit (Angled)	(29)
MPF Composite Application Tools Set 1	(30)
MPF Composite Application Tools Set 2	(31)
Replacement Tips	(34-35)

03 PG 36-45**Palettes**

MPF Composite Warming Palette	(39)
Marble Ceramic Palette	(40)
Glass Ceramic Palette	(41)
CAD Artistry Palette	(43)
Stain/Opaque Palette	(44)
Evolution Master Ceramic Palette	(45)

04 PG 46-49**Carving Instruments**

MPF Revolution Ceramic Sculpturing Instruments	(48)
Optimum™ Ceramic Carving Instruments	(49)

05 PG 50-71**Tools & Accessories**

MPF Zirco Texturizer	(52)
MPF Zir Grip	(52)
MPF Ceramo Blade	(53)
MPF Zirco Cleaner	(53)
MPF Zirco Scorer	(54)
Scribe™ Instrument by Joshua Polansky	(54)
MPF Ceramomix™ Zirconia Spatula	(55)
MPF Agate Mixing Spatula	(55)
MPF Debubblizer	(57)
MPF Master Klip	(58)
MPF Rainbow Drops Ceramic Colorants	(58)
MPF Brush Hydrator	(59)
MPF Brush Pool	(59)
Platinum Foil Technique by Ricardo Schäfer	(60)
MPF Drill Stands	(61)
MPF Surgical Drill Stands	(63)
MPF Implant Drivers	(65)
MPF Honeycombs & Pins	(66)
MPF Veneer & Crown Stabilizing Paste	(67)
MPF Peg Fix	(67)
MPF Waxer Instruments	(68)
MPF Crown Holders	(69)
MPF Spin Fix	(70)
Optimum™ Instrument Stand	(70)
MPF Texture Eyes	(70)
MPF All-Bright Diamond Paste	(70)
MPF Cube Stands	(71)
Liquid Dispencers	(71)
MPF Ceramic Table Mat	(71)

06 PG 84-91**Dental Training School**

Dental School Tools	(86-89)
Dental School Kit	(91)

01


Hybrid Series




MPF Hybrid Brush Tip

Patent pending (A.N. : CY202200003)


Due to the decrease of High-Quality Kolinsky, our solution was to incorporate a Felt Nib insert within the hair bundle. This absorbent Felt Nib when wet withholds water and due to the fact that the hair bundle is wrapped around the felt nib creates a constant humid environment. This keeps the hair in a constant humid state.

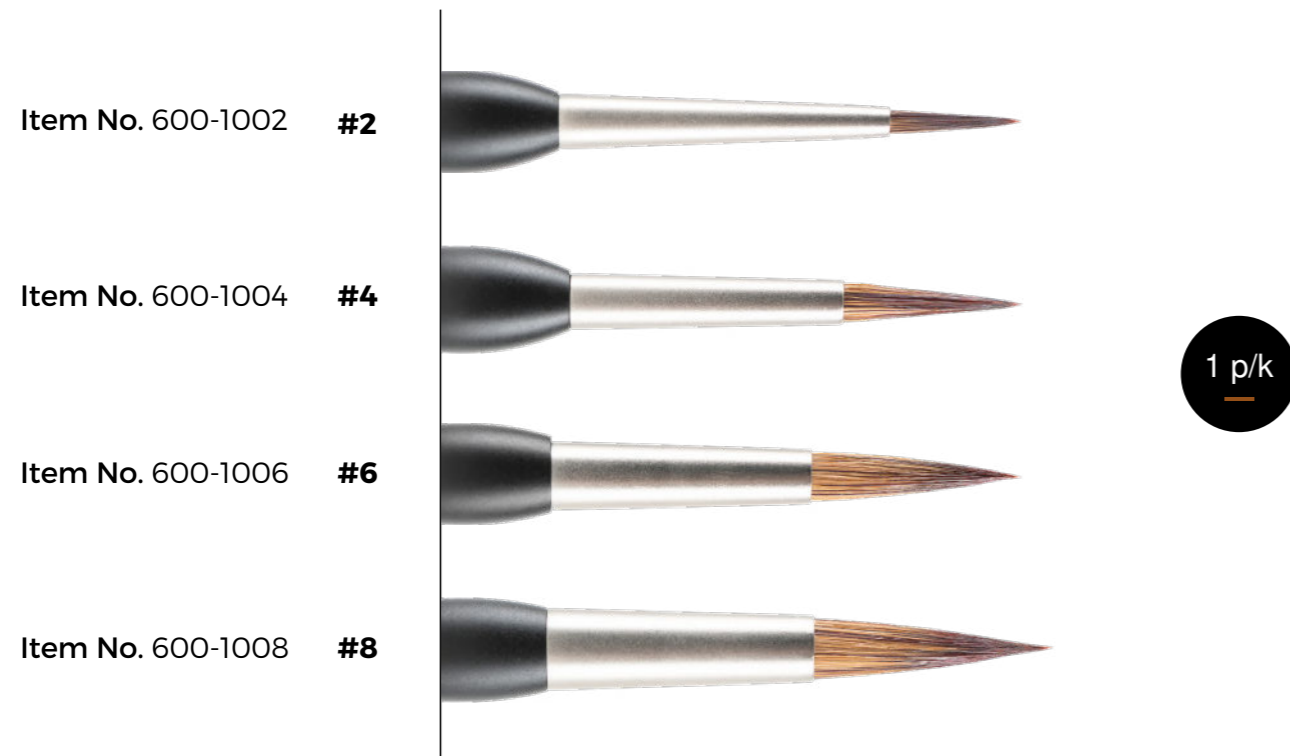

Absorbent Felt Nib

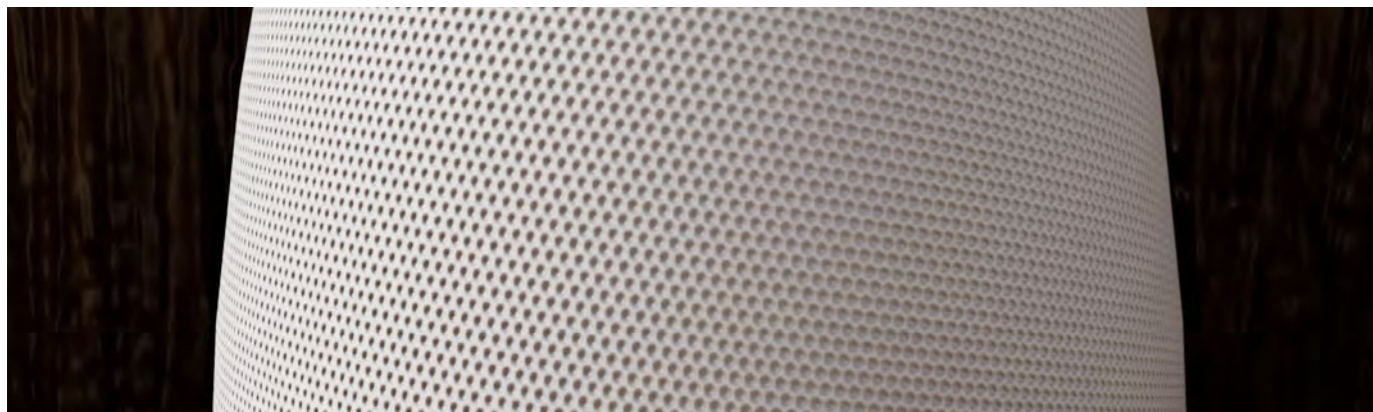

Compatible with every MPF Brush Handle.

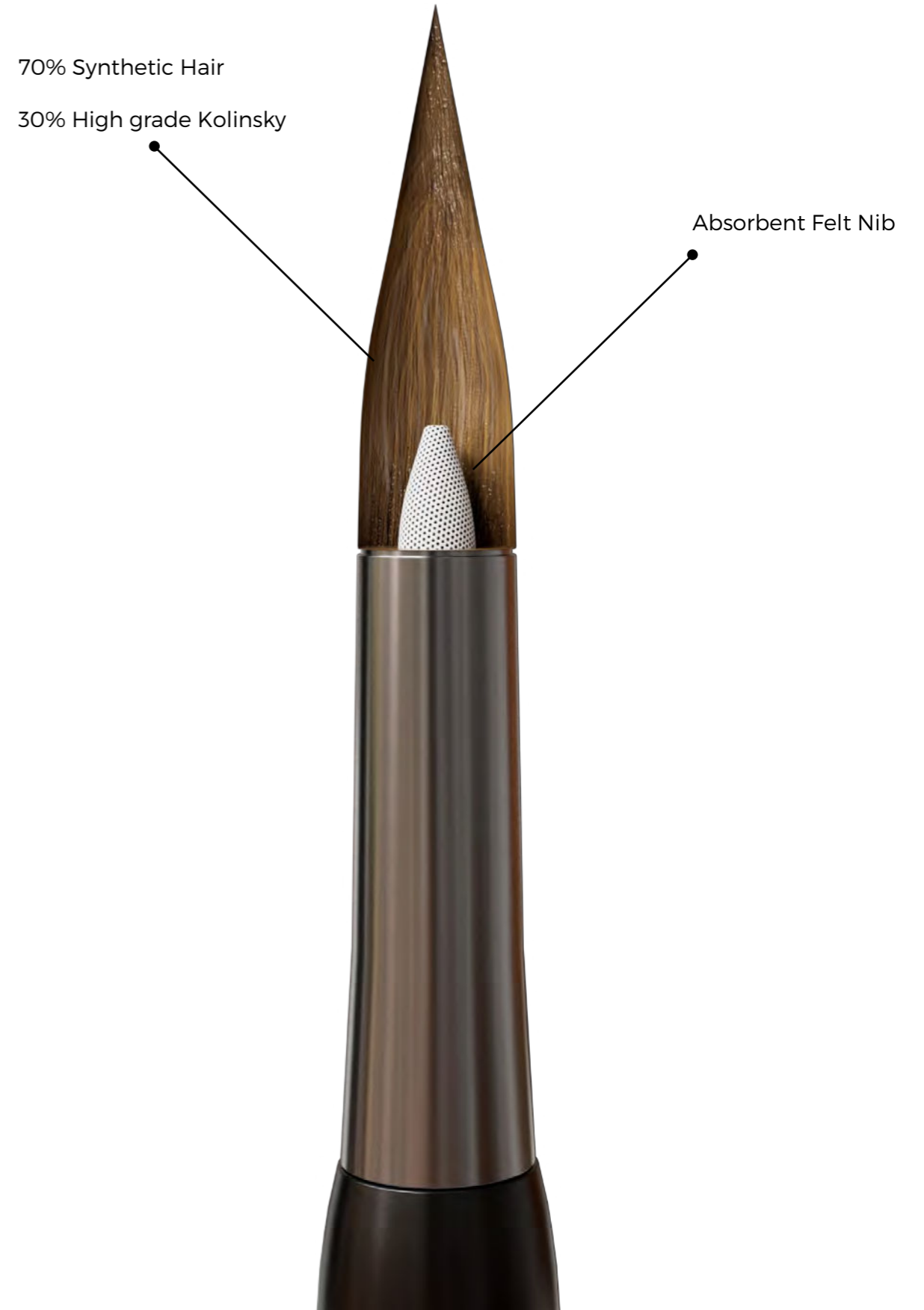


MPF Hybrid Plus Brushes & Brush Kit

Our recent Innovation in the dental brush industry, the Hybrid brush line with its absorbent internal nib implementation within the brush head, allowed us to create a much better humid environment for Synthetic hair.

The Hybrid Plus brush series are made out of 100% high grade synthetic hair that is perfect when being used for dental ceramics and staining. The absorbent material in the core of the hair bundle mimics the action of the natural hair, and delivers the exact humidity levels that are needed.


This innovative brush tip design is housed in our complete plastic handles, that also contain a few metallic balls within the ferrule of the brush tip. These balls are allowing more condensation when activated with movement.

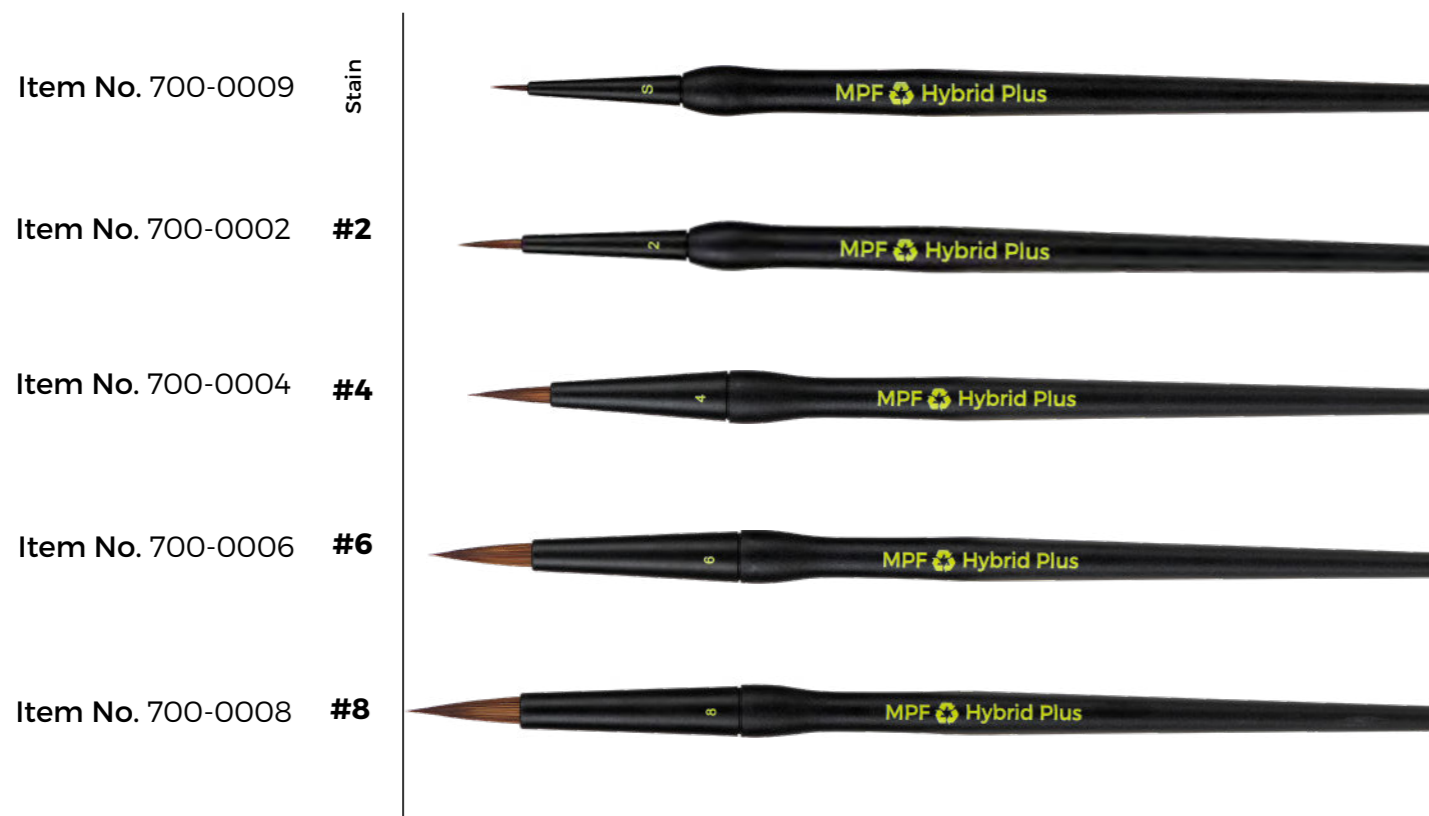

Item No. 700-0000
MPF Hybrid Plus 4 in 1 Brush Kit

100% Synthetic

Absorbent Felt Nib

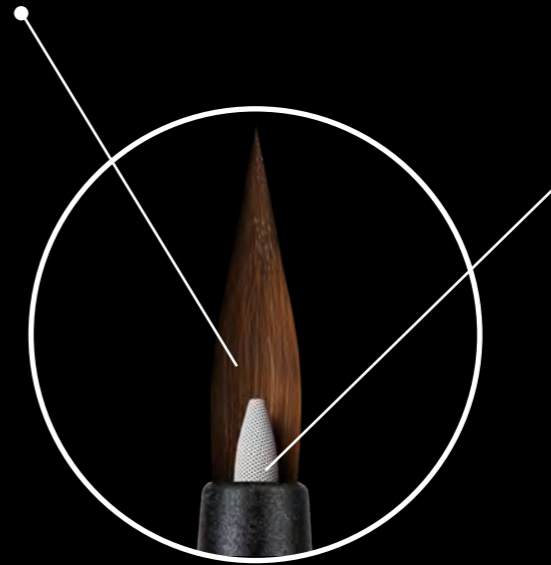

02

Brushes & Brush Kits



Optimum™ Spring Ceramic Brushes & Kits: Mens

Spring activated tip. Perfectly balanced. Replaceable premium Kolinsky sable brush tips. Integrated thorn allows to release pinned die-segments from model base.



Item No. 102-2000

Optimum™ Spring 3 in 1 Brush Kit.
Citrus

Item No. 102-3000

Optimum™ Spring 3 in 1 Brush Kit.
Champagne

Item No. 102-4000

Optimum™ Spring 3 in 1 Brush Kit.
Midnight Blue

Optimum™ Spring Ceramic Brushes & Kits: Ladies

Optimum™ Spring Ceramic Brushes Kits contain all 3 Brush Tip sizes (#4, #6, #8), plus the handle of your choice.



Item No. 103-1000

Optimum™ Spring 3 in 1 Brush Kit.
Pink

Item No. 103-2000

Optimum™ Spring 3 in 1 Brush Kit.
Purple

Item No. 103-3000

Optimum™ Spring 3 in 1 Brush Kit.
Midnight Blue

Optimum™ Master Brushes & Kits

All the tools you need in one package. Replaceable premium Kolinsky sable brush tips. Integrated thorn allows to release pinned die-segments from model base.

The Kit Contains: Stain, Opaque, Contour brushes, plus a Builder Brush of your choice.

Item No. 102-0009

Stain



Item No. 102-0010

Opaque



Item No. 102-0011

Contour



Item No. 102-5004

Optimum™ Spring Master Kit.
Size 4



Item No. 102-5006

Optimum™ Spring Master Kit.
Size 6



Item No. 102-5008

Optimum™ Spring Master Kit.
Size 8

* For Replacement Tips, see page 32

MPF Revolution Brushes & Kits

Light-weight, ergonomic aluminium handle. Replaceable premium Kolinsky sable brushtips.

Revolution Same 3 tips Brush Kit contains: 3 Rplc Tips of the same size.
Available size Tips: #2, #4, #6, #8



Item No. 200-3002
MPF Revolution Same 3 tips
Brush Kit. **Midnight Black**
Size 2

Item No. 200-3004
MPF Revolution Same 3 tips
Brush Kit. **Silver Moon**
Size 4

Item No. 200-3006
MPF Revolution Same 3 tips
Brush Kit. **Lead Grey**
Size 6

Item No. 200-3008
MPF Revolution Same 3 tips
Brush Kit. **Sea Blue**
Size 8

MPF Revolution Master Brushes & Kits

Revolution Master Brush Kit contains: Stain, Glaze, Contour, brushes plus a Builder Brush.

Available size Tips: #2, #4, #6, #8



Item No. 200-2002
MPF Revolution Master
Brush Kit.
Size 2

Item No. 200-2004
MPF Revolution Master
Brush Kit.
Size 4

Item No. 200-2006
MPF Revolution Master
Brush Kit.
Size 6

Item No. 200-2008
MPF Revolution Master
Brush Kit.
Size 8

MPF Revolution Handles


Light-weight, ergonomic aluminium handle.

Compatible with every MPF Replacement Brush Tip.

Item No. 200-1



Item No. 200-2

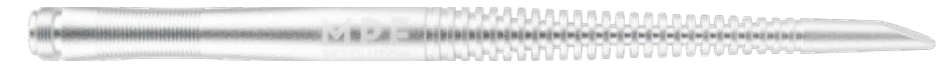

Item No. 200-3



Item No. 200-4



Tail Detail



Jungo Endo Signature Brushes

What makes this brush stand out from the others is its long premium kolinsky brush tip. Stay close to nature series was inspired by Jungo Endo a Japanese dental technician located in USA.

Item No. 118-0000



Item No. 118-1000



Item No. 118-1500



MPF Signature Brush

Item No. 118-2000



Stipple by Stefan Picha

A stippling brush for creating effects on the gingival surfaces.

Item No. 118-3000



MPF Standard Brushes & Kits

Innovative spring activation system. The perfect tip with a simple 'flick'. Replaceable premium Kolinsky sable brush tips.

All Standard 3 in 1 Brush Kits contain: all 3 brush tip sizes.
Available sizes #4, #6, #8

Item No. 100-0004 **#4**



Item No. 100-0006 **#6**



Item No. 100-0008 **#8**



Item No. 100-2000
MPF Standard 3 in 1 Brush Kit.
Silver



Item No. 100-3000
MPF Standard 3 in 1 Brush Kit.
Blue



Item No. 100-4000
MPF Standard 3 in 1 Brush Kit.
Black

MPF Zirconia Stain Brushes & Kit

A revolutionary solution for your stain needs. Metal and rust free. Plastic handle and ferrule. Synthetic hair. Ideal for water based liquids and acids. Non-replaceable brush tips.

Item No. 105-1001 **#1**



Item No. 105-1002 **#2**



Item No. 105-1003 **#3**



Item No. 105-1000
MPF Zirconia Stain 3 in 1 Brush Kit.

MPF Synthesis Brushes & Kit

3 p/k

The Synthesis Series Philosophy is to create a more humid environment to the synthetic hair by adapting Kolinsky into its core.

The humid properties that the kolinsky core encases are distributed through the longer synthetic hair, creating a wonderful combination at a very affordable price.

Item No. 105-2000 **#2**



Item No. 105-4000 **#4**



Item No. 105-6000 **#6**



Item No. 105-8000 **#8**



Item No. 105-0000
MPF Synthesis 4 in 1 Brush Kit



MPF Classic Kolinsky Brushes & Kit

Premium Kolinsky sable brush tips. New plastic handle and ferrule. Tip de-humidifying system (internal ball bearings.)

MPF Classic Kolinsky Brush Kit contains sizes: #2, #4, #6, #8

Item No. 105-0002	#2	
Item No. 105-0004	#4	
Item No. 105-0006	#6	
Item No. 105-0008	#8	
Item No. 105-0009	Stain	
Item No. 105-0010	Opaque	
Item No. 105-0011	Contour	
Item No. 105-0012	Glaze	




Item No. 105-0001
MPF Classic 4-1 Brush Kit

MPF Composite Stain Brushes STRAIGHT

The MPF COMPOSITE BRUSH series have been created with a special fiber optic filament that is very fine in texture but also firm to assist in application.



MPF Composite Stain Brushes ANGLED

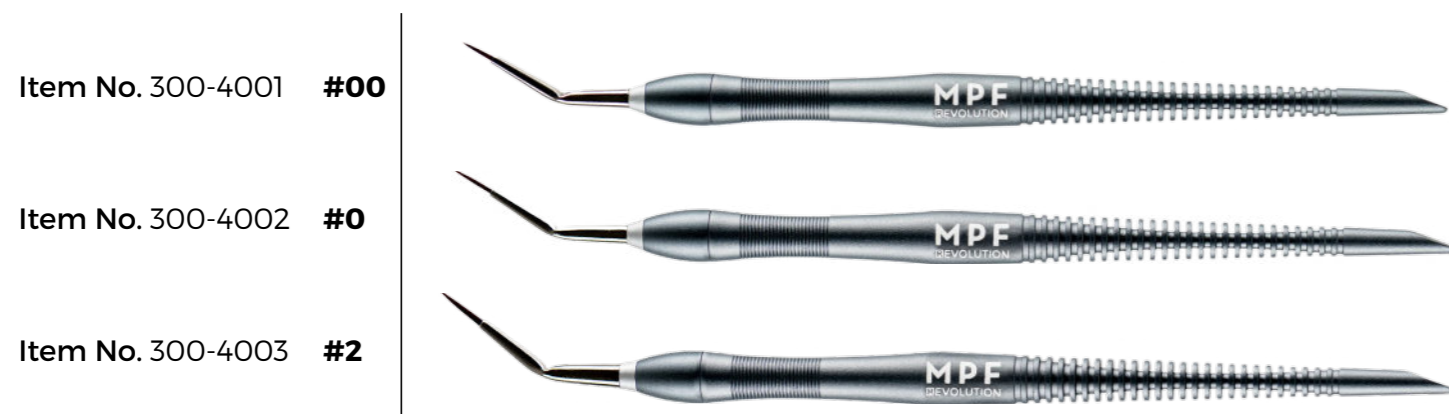

Item No. 300-3000
MPF Composite 3 in 1 Stain Brush Kit.
Straight

Item No. 300-4000
MPF Composite 3 in 1 Stain Brush Kit.
Angled

MPF Composite Flat Brushes STRAIGHT



MPF Composite Flat Brushes ANGLED




Item No. 300-5000
MPF Composite 3 in 1 Flat Brush Kit.
Straight

Item No. 300-6000
MPF Composite 3 in 1 Flat Brush Kit.
Angled

MPF Composite Application Tools (Direct) Set 1

MPF Direct composite applicators are designed with ergonomics, comfort and balance in mind. The unique design of each and every tip allow for intraoral application and also for the dental technicians lab composite applications as well.

The handles are made from anodised aluminium and the tips are stainless steel that have a special BLACK DIAMOND coating.

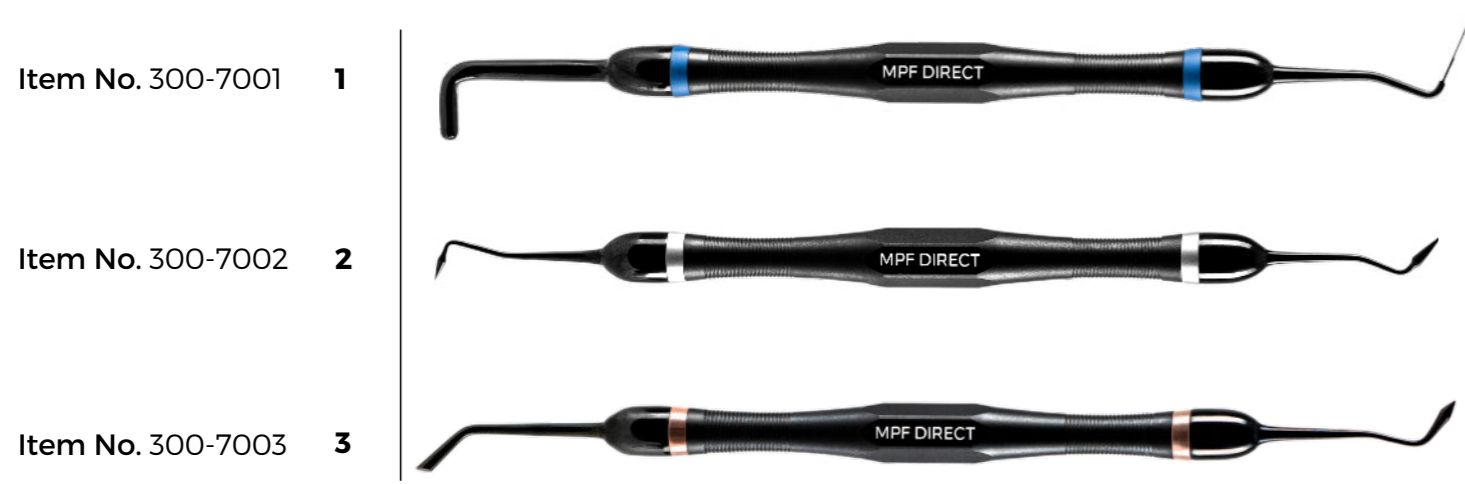

Item No. 300-7000
MPF Composite Application Tools (Direct)
Set 1

MPF Composite Application Tools (Direct) Set 2


These tools are all autoclavable and can withstand any sterilisation process.

They are NOT to be cleaned by using wire brush techniques as this can destroy the coating process. Care must be taken by using a soft brush to remove any reminiscence of composite material.



Item No. 300-8000
MPF Composite Application Tools (Direct)
Set 2

Optimum™ Replacement Tips

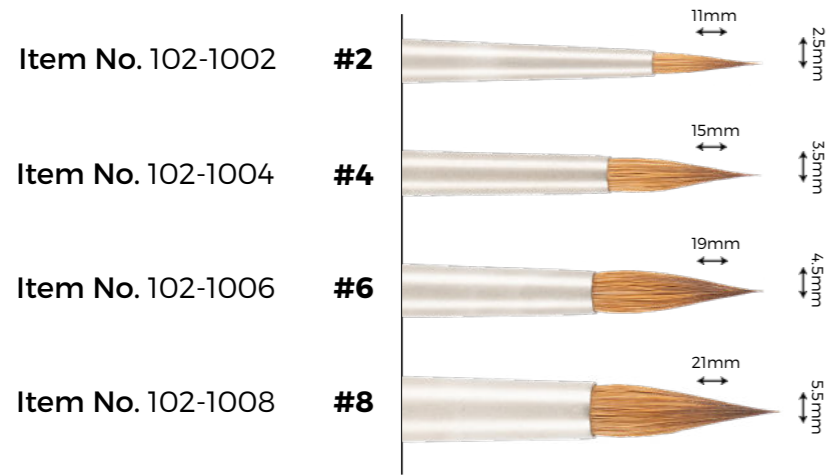


2 p/k




2 p/k

MPF Standard Replacement Tips

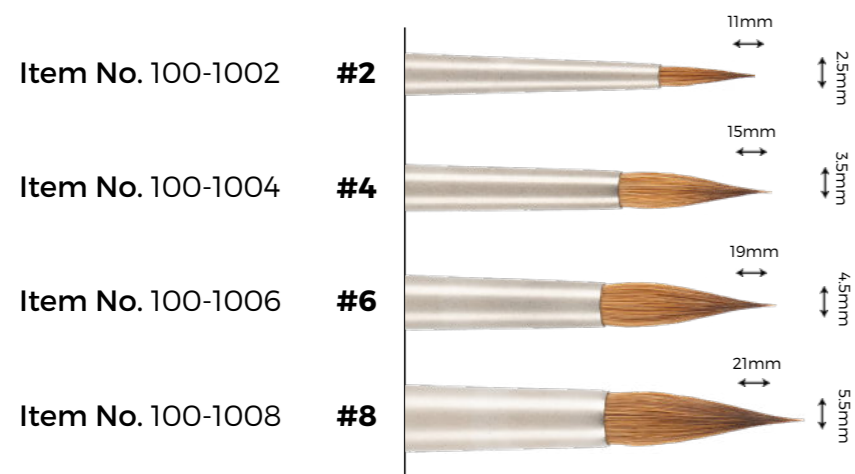


1 p/k

MPF Revolution Replacement Tips

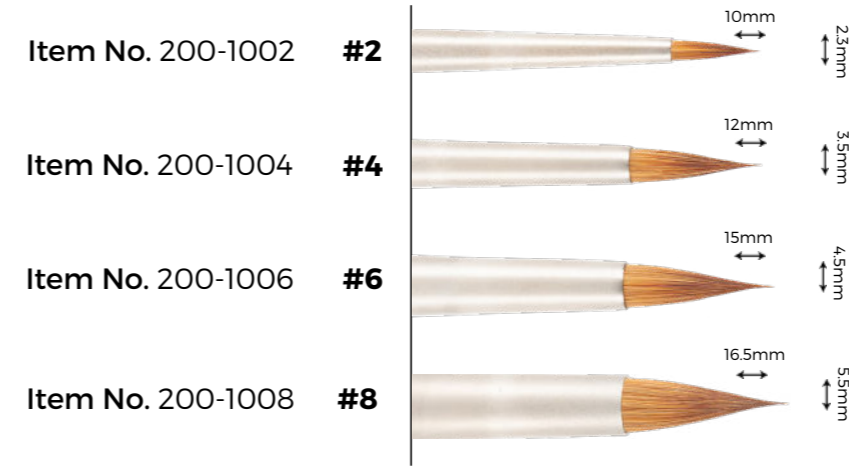


1 p/k

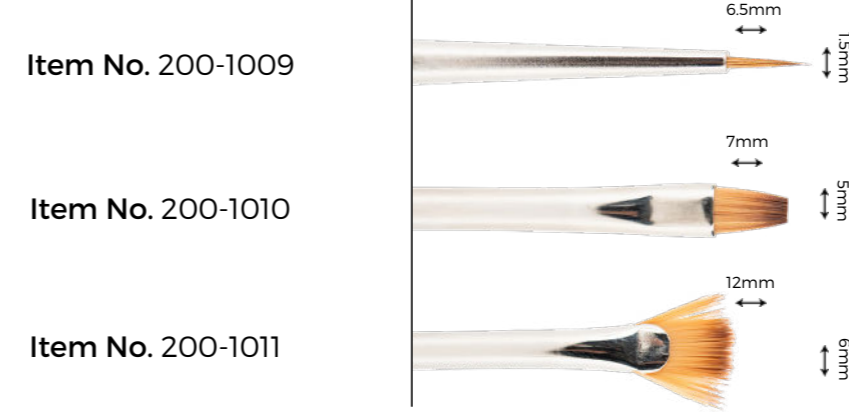


1 p/k

MPF Signature Replacement Tips

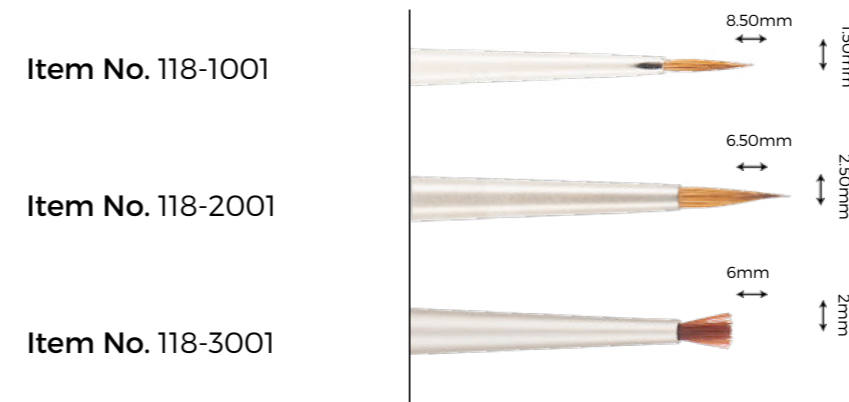


2 p/k

MPF Composite Stain (Straight) Replacement Tips

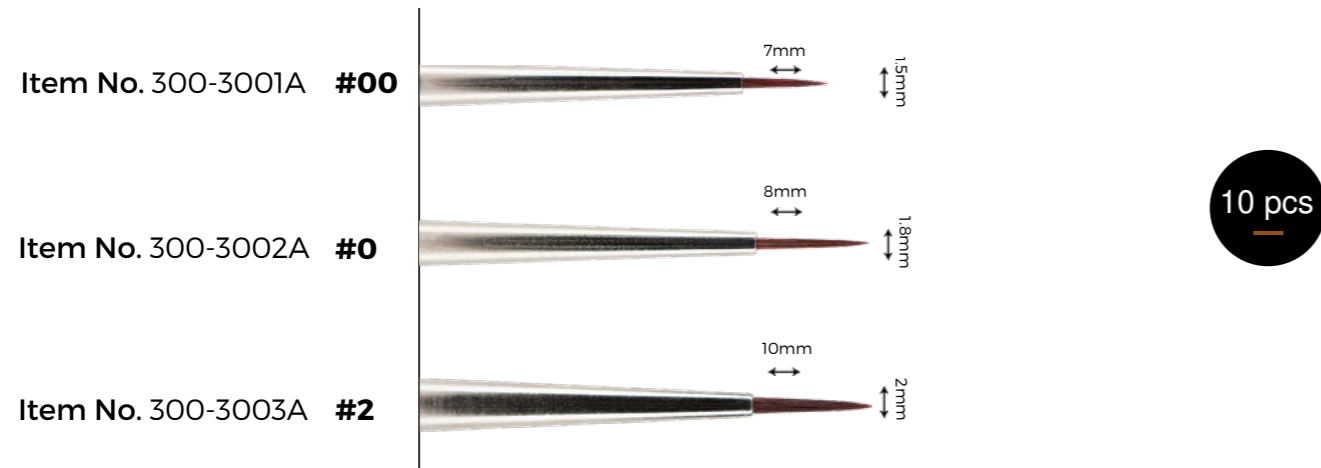


MPF Composite Flat (Straight) Replacement Tips

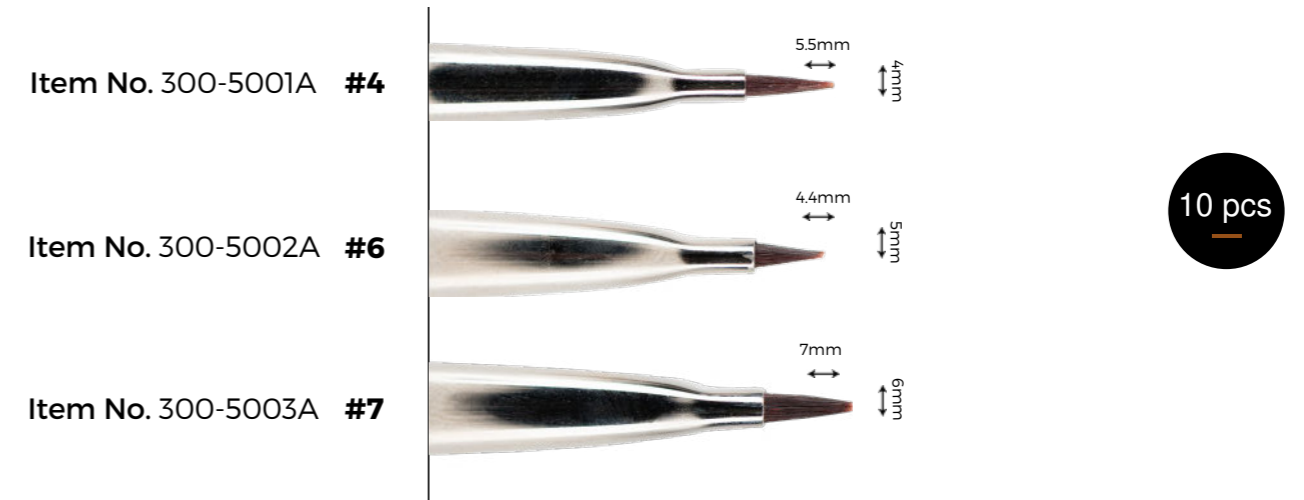


MPF Composite Stain (Angled) Replacement Tips



MPF Composite Flat (Angled) Replacement Tips

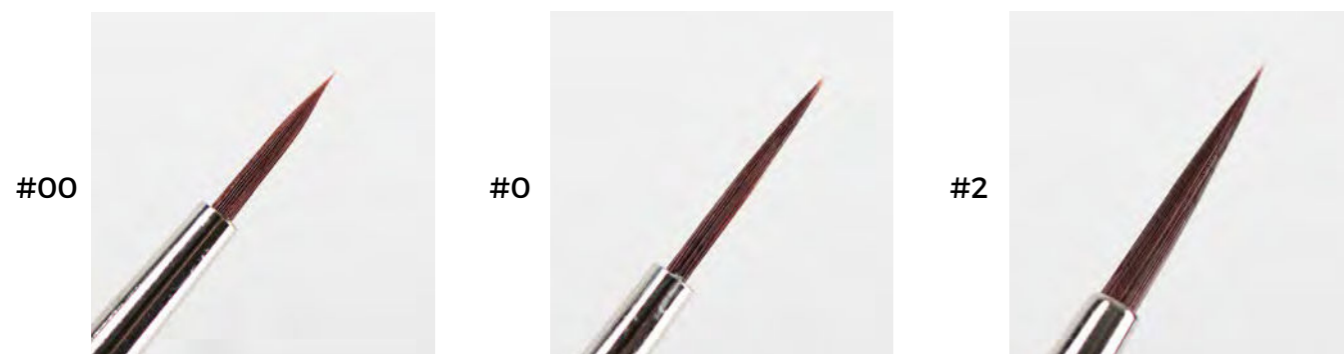

Close-up

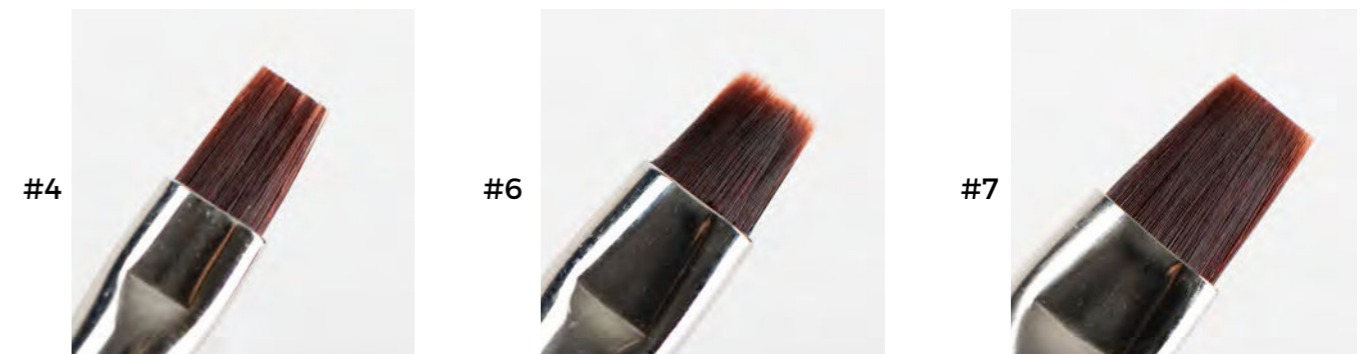

Close-up

03

Palettes



MPF Composite Warming Palette

Introducing the first ever MPF Composite Warming Palette. This unique innovation has been designed to bring the composite to a more malleable consistency which assists the user in manipulating the composite when restoring small or long span restorations especially when the need for extensive gingival adaptation in hybrid cases is required. The palette case has been created from a polycarbonate material that is called "Makrolon", and has stable UV properties that prevents direct UV light from polymerizing the composite within the tray. Furthermore, the lid of the palette box has the ability to remain open at any height you desire, in order to prevent direct light from your laboratory surroundings.

Item No. 115-9000

Includes: Power plug cable, blade scraper.

The glass plate is warmed via an element adhered to the bottom of the glass plate which is connected via a plastic connector that is connected to the electronics on the main body of the palette. The glass plate is 6mm thick, to avoid the risk of breaking at different temperatures.

The palette features an On/Off button on the top right of the inner palette. The same button can change the temperature LCD display from Celsius to Fahrenheit by simply clicking it. When holding down the On/Off for 3 seconds the palette turns Off.

The LCD display depicts the temperature of the actual glass plate either in Celsius or Fahrenheit and can be risen or lowered via the +/- buttons. When the desired temperature has been chosen, the palette will need a few minutes in order to reach the desired temperature levels.

The palette comes with a blade scraper which is used to clear the glass plate of composite debris without scratching the glass surface. Once the composite reminisce is removed from the glass with the scraper, the glass surface can be wiped down with a clean paper towel and alcohol to clean the surface.



The palette will automatically turn off after 5 hours in case the unit has been left on by mistake.

Marble Ceramic Palette

Italian marble surface with nylon worksheet.
Durable base. Capillary action keeps ceramic mix continuously wet.

Item No. 115-1000

Includes: marble tray, 2 worksheets, liquid dispenser

The unique marble used to create this palette allows controlled amounts of water to permeate and wet the nylon worksheet, ensuring a uniform, consistent and workable ceramic mix. The nylon worksheet covering the marble is kinder to Kolinsky brushes and helps keep a finer point for longer periods of time.




Glass Ceramic Palette

Glass work surface. Durable base. Synthetic leather strips keep ceramic mix continuously wet.


Item No. 115-0000

Includes: glass tray, spatula, 4 strips, liquid dispenser

The Glass Ceramic Palette, created for ceramists who prefer a glass working surface, uses a system of synthetic leather strips to keep ceramic mix continuously wet. Water from the base permeates up through the strips and into the ceramic stack, creating a homogeneous ceramic mix.

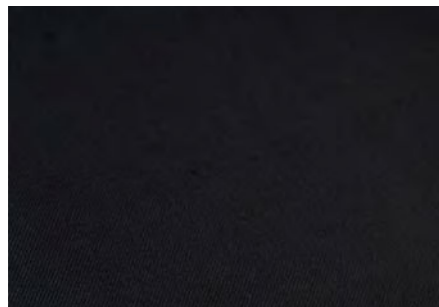

Item No. 115-1001
Nylon Palette Worksheets 2/pk

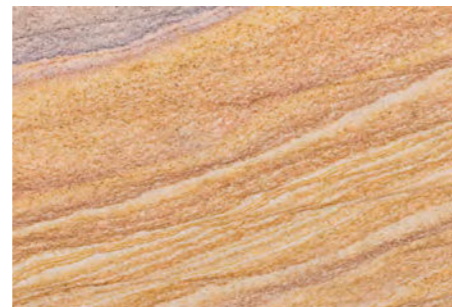

Item No. 115-1002
Marble Stone



Item No. 115-0001
Palette Absorption Strips 4/pk



Item No. 115-0002
Replacement Glass

CAD Artistry Palette



Moisture and consistency control for ceramic material. Dual working platforms. CAD-CAM Technology. This palette is the perfect solution for ceramists working with CEREC/IN-LAB and CAD-CAM technology.

Item No. 115-7000

Includes: tray, 2 strips, 1 worksheet, sponge, liquid dispenser

The CAD Artistry Palette was designed with two different options of humidification to keep the proper moisture and proper consistency of ceramic material. The two platforms on the ceramic tray allow you to choose between humidification with water or with liquid.

Choose between working on a nylon sheet or using the synthetic strips. Eleven stain compartments act as storage for your favourite stains and colours.

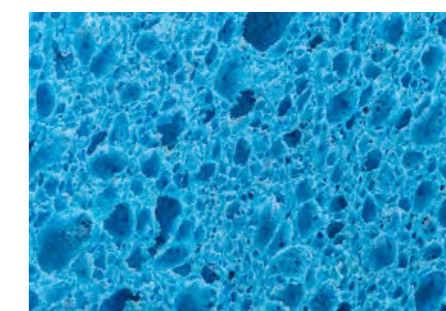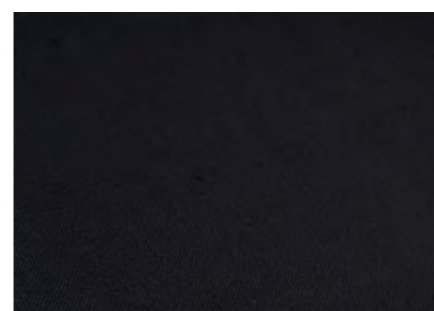

Item No. 115-7001
Nylon Palette Worksheet

Item No. 115-0001
Palette Absorption Strips 4/pk

Item No. 115-4005
Replacement Sponge

Stain/Opaque Palettes

Smart storage solution for stain and paste opaque materials. Durable plastic structure. Tight seal lid.

Item No. 115-5000
Black

Item No. 115-6000
White

The unique, angled, step-type design of the Stain & Paste Opaque Palette allows easy access to materials. Reservoirs have a very unique design, elongated towards the rear so that material flows to the back of the cup, preventing spillage. The tight seal lid prevents oxygen drying out materials. A colour chart helps keep material colours in order.



Evolution Master Ceramic Palettes

Glass work surface. Synthetic leather strips keep ceramic mix continuously wet. Concealed stain palette draw.

Item No. 115-4000
Black

Item No. 115-4100
White

Includes: glass tray, 4 strips, 1 nylon worksheet, sponge, stain brush

The Evolution Master Ceramic Palette has a sleek, contoured design. It's ergonomic design is user friendly, allowing easy access to your ceramic materials without obtrusive, robust parts inhibiting your movements to the anterior portion of the tray.

Item No. 115-4006
Replacement Glass Tray (Black)



Item No. 114-4007
Replacement Glass Tray (White)



Item No. 115-4001
Replacement Nylon Worksheet



Item No. 115-4003
Replacement Synthetic Strips 4/pk



Item No. 115-4005
Replacement Sponge




04

Carving Instruments




MPF Revolution Ceramic Sculpturing Instruments

Precision ceramic carving tools. Perfectly balanced. Replaceable tips. Ergonomic handles.

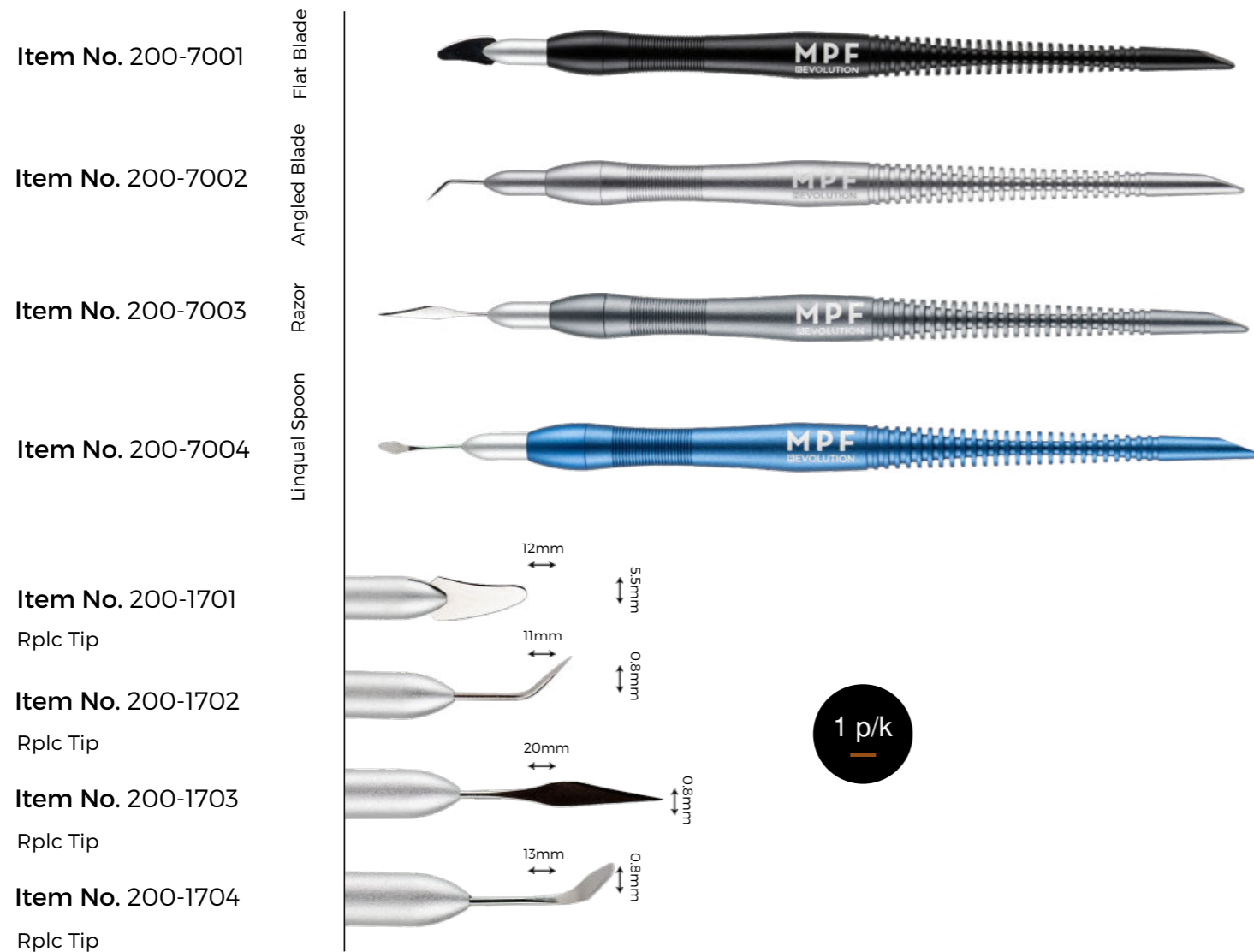


Optimum™ Ceramic Carving Instrument

Integrated thorn allows to release pinned die-segments from model base.

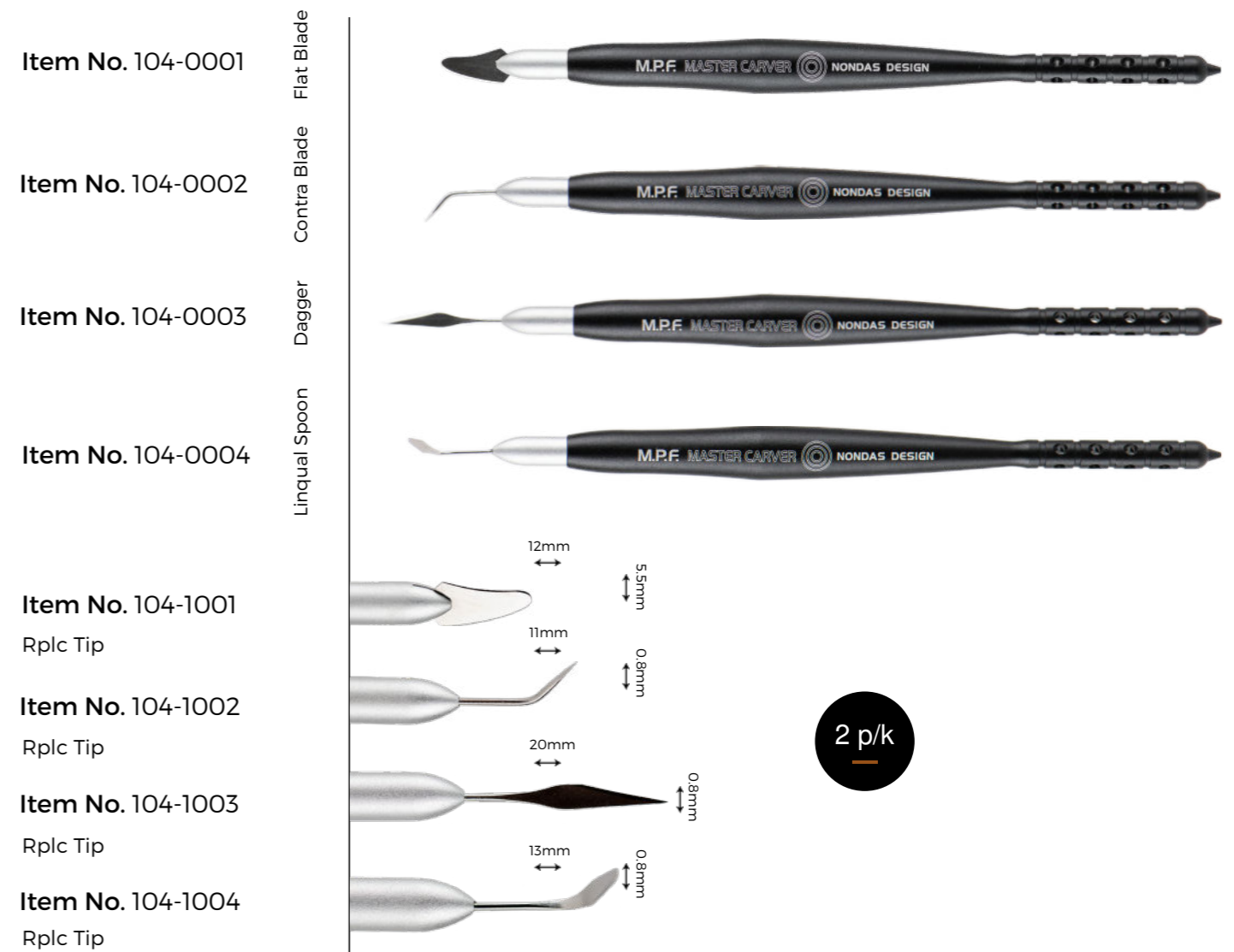

05

Tools & Accessories

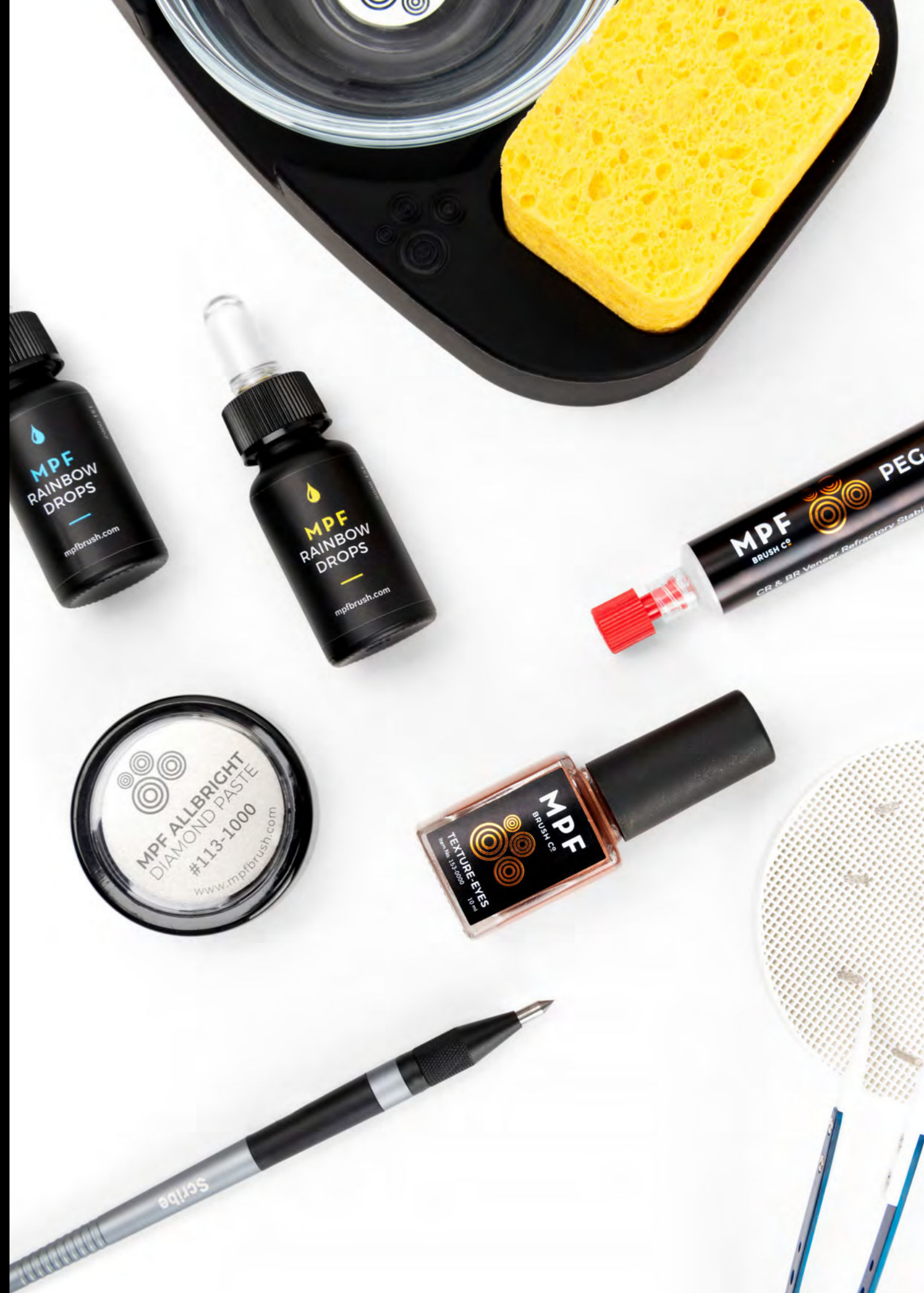

MPF Zirco Texturizer

MPF Zirco Texturizer allows you to create soft subtle textures on unsintered Zirconia.

Item No. 400-1000

Before the sintering process this stainless steel wire brush can be used to create asymmetric lines that look very natural after the sintering process, and also remove any fine zirconia powder that may be left on the outer surface of the restoration.



MPF Zir Grip

MPF Zir Grip Tweezers have been designed to work with all kinds of ceramic & zirconia restorations in mind.

Item No. 400-5000

These ergonomically designed tweezers boast tips made from sintered Zirconia and are shaped in such a way to assist the ceramist when positioning his crowns, veneers and small bridges onto the refractory pins before firing. These tips being from zirconia avoid leaving any dark marks that the current metallic tweezers create.



MPF Ceramo Blade

The MPF WCeramo Blade can work directly on unsintered Zirconia creating intricate surface desing.

Item No. 400-3000

Item No. 400-3001
Replacement Tip

This MPF Ceramo Blade can work directly on unsintered Zirconia creating intricate surface desing such as crack lines & deep fossa in the occlusal region, the superior sharpness allows for the ease of removing portions of zirconia to create wear facets & creative details before your final sintering.



MPF Zirco Cleaner

This newly designed brush tip has been created to allow for ease of cleaning within zirconia milled crowns, bridges & implant restoraions.

Item No. 400-4000

Item No. 400-4001
Replacement Tip

Zirconia milling creates a very fine powder after milling, and is difficult to not only see, but to remove, so this new brush tip has been designed with very stable synthetic hair in a very tightly shaped bundle, which gives excellent properties for removing this fine powder from your work. the fitting surface.



MPF Zirco Scorer

MPF Zirco Scorer allows you to create amazing surface texture on unsintered zirconia.

Item No. 400-2000

Item No. 400-2001
Replacement Tip

Due to the fact that Zirconia is a sensitive material Zirco Scorer allows us to create amazing surface texture on unsintered zirconia, pit fissure grooves and crack lines are all possible with this amazingly sharp scoring instrument and being made of Zirconia prevents any contaminations to the raw material.



MPF Ceramomix™ Zirconia Spatula

Ceramic mixing instrument that prevents contamination of ceramic material.

Item No. 108-1000

The spatula has been developed and combined with our spring activated technology, to relieve pressure on the zirconia blade when mixing, this has enhanced the life span of this product by 10 fold.



Scribe™ Instrument by Joshua Polansky

Amazing surface texture and crack lines are all possible with this amazingly sharp scoring instrument.

Item No. 107-0000

Item No. 107-0001
Replacement Tip



MPF Agate Mixing Spatula

Made from Brazilian Agate Stone, with the tail end shaped at an angle for easy extraction of ceramic powder.

Item No. 108-0000



MPF Debubbler

The MPF Debubbler is an advanced formula, clinically designed to reduce surface tension, thus eliminating bubble formation when pouring dental model casts from impressions and investing wax patterns. A simple mist spray into silicone, rubber base or alginate impression materials relieves surface tension that our impression materials create and this leads to a flawless replication of the impression in the thinnest and most detailed of areas. This result applies to both model replications and also wax castings.

Item No. 157-0000
300ml

Item No. 157-1000
1000ml

The MPF Debubbler comes in a 300ml bottle and its Spray Nozzle and also 1000ml refill bottles.



Key Benefits

Reduces surface tension on wax and impressions.

Allows gypsum and investment to flow smoothly into all of your impressions.

Ensures bubble free models and smooth precision castings.

Guarantees consistent, accurate results.

MPF Master Klip

The Master Klip Tweezers are a spring activated hand tool that has stainless steel pointed holders that are activated and deactivated when the base of the tweezer is pressed.

Item No. 112-1000

These Tweezers are very useful for holding onto delicate prosthesis such as crown and veneers when staining or creating additions.



MPF Rainbow Drops Ceramic Colorants

High quality colourants for your ceramic mix .

Item No. 151-0000 (Kit)

Item No. 151-0001 (Red)

Item No. 151-0002 (Blue)

Item No. 151-0003 (Orange)

Item No. 151-0004 (Brown)

Item No. 151-0005 (Yellow)

Item No. 151-0006 (Violet)



MPF Brush Hydrator

MPF Brush Hydrator was made to keep you brushes perfectly hydrated throughout the day.

Item No. 109-0000

Pour some water in the glass reservoir, and simply mount the brush on the hydrator for about 5 minutes.

After 5 minutes of soaking, the water permeates its way into the hair bundle, giving perfect moisture level required for application.

ATTENTION: Hydrator is not a storage solution.



* brush not included

MPF Brush Pool

Ergonomic and durable design for cleaning and storing your favourite MPF brushes. The base is made out of silicon which makes it very easy to clean and dry.

Item No. 111-0000
Brush Pool

Item No. 111-0001
Replacement Bowl

Item No. 111-0002
Replacement Sponge



Platinum Foil Technique Kit by Ricardo Schäfer

In close collaboration with one of our okey opinion leader, Ricardo Schäfer from Argentina, we developed an adapter tool that boasts a double ended tool made from aluminium that has two Zirconia adaptation tips, designed by Ricardo Schäfer

Item No. 117-0000



The Kit contains:

- Adapter Tool: these ergonomic tips have been shaped to assist in the forming of the platinum foil onto the die.
- Scissors: a pair of sharp scissors that are designed to cut the platinum foil with ease.
- Master Klip: a tool used for holding onto your dies when staining and glazing.
- Tweezers: tweezers that are designed with non invasive tips to assist in the removal of the platinum foil after the final glaze firing.



MPF Drill Stands

The MPF Drill Stand allow us to store our burs according to the task at hand. The aluminum base is quite heavy and has non-slip silicone feet for stability. Universal fit for standard Rotary Instruments & F.G Burs.



Item No. 155-0001



Item No. 155-0002



Item No. 155-0003



Item No. 155-0004

MPF Surgical Drill Stands


Our drill stand is made from Aluminum and is color anodized in three radiant exciting colors. The three level design allows for a much easier visual of the head of the drills due to the fact that they are on different levels.

On the level closest to the base you can position up to 26 friction fit burs, on the middle level you can position up to 19 friction fit burs, and on the upper level you can position up to 12 regular hand motor burs.

The Surgical drill stand is connected to a flat base that allows the top half to spin freely for quick and easy access when choosing the appropriate bur for the task at hand.

The MPF Surgical Drill stand is fully autoclave-able.

The Drill Stand should not be put into an ultrasonic with chemicals !

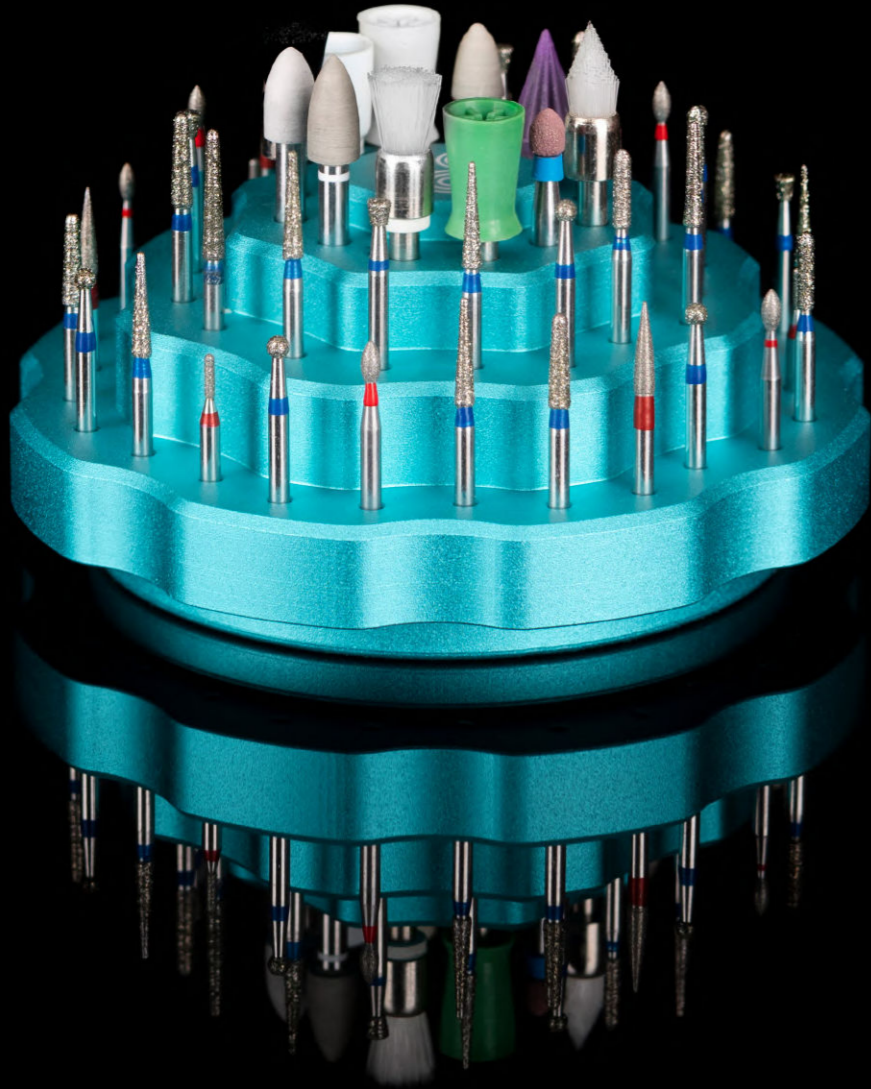

Item No. 156-0000

Aqua

Item No. 156-0001

Citrus

Item No. 156-0002

Rosegold



*Drills are not included

MPF Implant Drivers

The screwdriver bits are of the highest quality tool steel and heat treated for maximum hardness and resilience to wear and tear. The ergonomic handle is made from aluminum and color anodized and numbered in 8 brilliant colors for easy identification.



Item No. 154-0001
Hexagon 0.9



Item No. 154-0002
Hexagon 0.95



Item No. 154-0003
Hexagon 1.2



Item No. 154-0004
Hexagon 1.27



Item No. 154-0005
Hexalobular T5.5



Item No. 154-0006
Hexalobular T6



Item No. 154-0007
Hexalobular Ball



Item No. 154-0008
Slot



Item No. 154-0000
Kit

Kit comes in an aluminum base for instant access and storage.

MPF Honeycombs & Pins

All new 'hot zone' technology. Zirconia pins with an elevated base for crowns and implants.

Item No. 114-1000
Honeycomb Tray with 10 Ceramic Pins

Item No. 114-1001
Ceramic Pins, 10 pcs

Item No. 114-2000
Honeycomb Tray 10 Metal Pins

Item No. 114-2001
Stainless Steel Pins, 10 pcs

Item No. 114-3000
Honeycomb Tray for Implant bound CR & BR with 4 ceramic pins

Item No. 114-3001
Custom firing Ceramic Pins for Implant bound CR & CB, 4 pcs.

Item No. 114-4000
Honeycomb Tray for Veneers 40mm x 20mm

Item No. 114-4001
MPF Alumina Veneer Pins 6pcs



MPF Veneer & Crown Stabilizing Paste

A stabilizing paste used to safely adhere metal and metal free prosthesis to your model dies with superior holding properties. Burns residue free.

Item No. 150-0001



MPF Peg Fix

A new, revolutionary, environmentally friendly refractory stabilizing paste made from PCW fibers, suitable for firing with metal or metal free prosthesis. Can be completely removed from your prosthesis with a short blast of compressed air, leaving it totally residue free.

Item No. 150-0002
(1 tube)

Item No. 150-0003
(3 tubes)



MPF Waxer Instruments

Lightweight, aluminium handle. Replaceable stainless steel blades. Two instruments in one.

Item No. 119-0001
Spoon/Fine Dropper.
Red



Item No. 119-0002
Large LeCron/Angled Spade.
Blue



Item No. 119-0003
Small LeCron/Large Dropper.
Purple



Item No. 119-1001
Spoon



Item No. 119-1002
Fine Dropper



Item No. 119-1003
Large LeCron



Item No. 119-1004
Angled Spade



Item No. 119-1005
Small LeCron



Item No. 119-1006
Large Dropper



Item No. 119-2000
Kit



MPF Crown Holders

Assistive tool to safely hold crowns in place while working, without chipping or marking. Special gum material is applied to the desired end (retention knurl or cut out), onto which the crown is positioned.

Item No. 152-0001
3-4mm Handle

Item No. 152-0002
2-2.5mm Handle



The MPF Crown Holder is an instrument designed to successfully hold crowns in place while working, without chipping or marking them.



Item No. 152-0003
White

Item No. 152-0004
Black



The adhesive gum is guaranteed not to leave sticky, greasy residue on the crown or on fingers.



Item No. 152-0000
MPF Crown Holder
Complete Kit



MPF Spin Fix



Item No. 106-0000

Ensures easy access and stability to crowns placed onto the tray, for ceramists to complete finishing touches with ease.

Optimum™ Instrument Stand



Item No. 104-4000

Keeps all your brushes and tools in one place. Easy Grab-N-Go design, durable and elegant.

MPF Cube Stands



Item No. 300-1000

Ergonomic stand for storing your brushes.
Blue

Item No. 300-2000

Ergonomic stand for storing your brushes.
Black

MPF Texture Eyes



Item No. 153-0000

An amazing copper paste that, when painted directly onto the surface of your crowns, veneers, bridge work or gypsum model, allows you to visualize and accurately replicate morphology and tooth structure.

MPF All-Bright Diamond Paste



Item No. 113-1000

Apply a small amount of paste after polishing ceramic or composite surfaces for a lustrous, high shine finish.

Liquid Dispensers



Item No. 115-0005
20 & 100 ml

Item No. 115-0003
Empty Modeling Liquid Dropper 100ml

Item No. 115-0004
Empty Stain Liquid Dropper 20ml

Ceramic Table Mat



Item No. 110-0000

06

Dental Training School Kit



MPF Dental Training School

The kit comes with 27 different hand tools that will be used without a doubt in the future career of every student. We have chosen high Quality tools that will stand the test of time.

Item No. 160-0001

Large & Small wax up knife for all waxing applications , with a knife shape and also a tapered scoop shape in the one handle.

MPF Lessman Knife Large



MPF Lessman Knife Small



Item No. 160-0002

Item No. 160-0003

Plaster knife with a large anterior blade for cutting and a leverage shape tip at the other end of the handle, usually used when deflasking Dentures from muffles.

MPF G.S Wax Knife



Item No. 160-0004

Straight edge plaster spatula for all mixing purposes.

MPF Plaster Spatula



Item No. 160-0005

Large & Small cement and acrylic mixing spatulas with a tapered tip and a flatter shaped rear tip, the handles have excellent grip due to the knurled surface of the handles.

MPF Cement Spatula Large



MPF Cement Spatula Small



Item No. 160-0006

MPF Dental Training School

Item No. 160-0007

Microblade holder with replacement blades used for removing bubbles from plaster models and general dental model preparation.

MPF Mini Blade Holder



Item No. 160-0008

Contra angle locking tweezers that can be used for holding sutures , carbon paper for bite registration purposes , general use in a Clinical or laboratory situation.

MPF Tweezer with Lock



Item No. 160-0009

Number 3 Scapel handle that is used with several blade sizes for a variety of applications in a laboratory or clinical situation.

MPF Scapel Handle



Item No. 160-0016

Straight scissors for the cutting of sutures, platinum foil and an abundance of applications in a clinical and laboratory situation.

MPF Scissors



Item No. 160-0017

Measuring calipers used for measuring the thickness of ceramic veneers and or metal framework, applied in both a clinical and laboratory situation.

MPF Decimal Caliper



MPF Dental Training School

Item No. 160-0010

Concave orthodontic pliers that can be used for several applications such as creating several types of Orthodontic clasps in a Laboratory or Clinical situation.

MPF Concave Plier




MPF Adam Plier




Item No. 160-0011

Adam pliers orthodontic pliers that can be used for several applications such as creating several types of Orthodontic clasps in a Laboratory or Clinical situation.

MPF Needle holder Straight

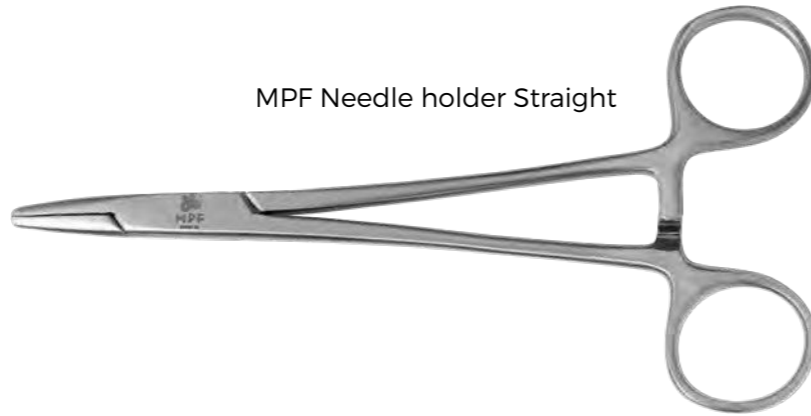

MPF Needle holder Curved

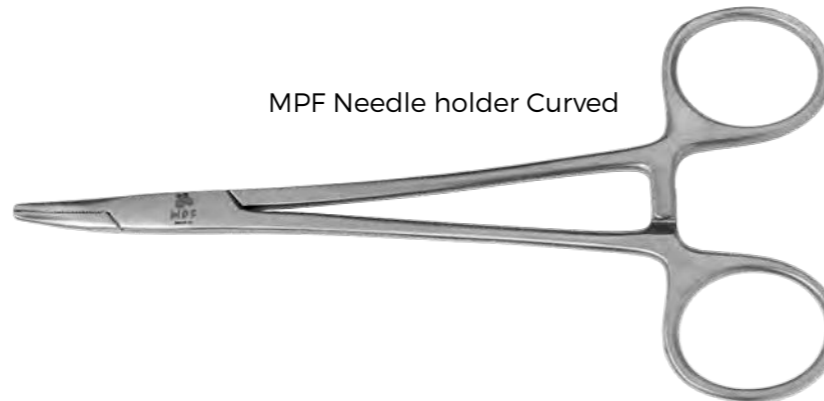

Item No. 160-0014

Needle holder Straight can be used in surgical procedures, but also used in a dental laboratory for several applications.

Item No. 160-0015

Needle holder Curved can be used in surgical procedures, but also used in a dental laboratory for several applications.

MPF Dental Training School

Item No. 160-0018

Alcohol burner which is used primarily in a clinical situation when working with wax or compound materials in a surgical situation.



Alcohol Burner

Item No. 160-0019

Alcohol torch can be used for the smoothing of wax surfaces in a clinical or laboratory situation.



Alcohol Torch

Item No. 160-0020

Silicon Plaster or Alginate mixing bowl which can be applied in both a clinical and laboratory situation.




Plaster Bowl

Item No. 160-0021

Silicon Bowls in small, medium and large sizes are primarily used for the mixing of powder to liquid materials which are applied in both a clinical and laboratory situation.

Silicone Bowls

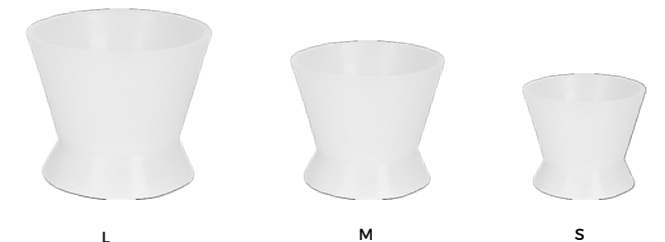

L

M

S

MPF Dental Training School Kit

Item No. 160-0000

Dental Training School Kit
*contains 27 tools in total.

Also included:

Combination Drill Kit

Ceramic Mixing Plate

MPF Classic Kolinsky Brushes #4 #6
(105-0004/105-0006)

Angle Spade/Large Lecron
(119-0002)

Large Dropper/ Small Lecron
(119-0003)

MPF Application Tools No2 & No4
(300-7002/300-8001)



Distributed in USA by:



WWW.ZUBLERUSA.COM
SALES@ZUBLERUSA.COM

Premium Dental Instruments

MPF BRUSH Co.

For more:



Distributed in USA by:



WWW.ZUBLERUSA.COM
SALES@ZUBLERUSA.COM

For more:

